



Parkell

DIGITEST SELF CONTAINED DIGITAL PULP TESTER




Item D626D
\$495.00

Specifically designed for patient comfort
Standard 9-volt batteries
Extremely economical
Simply place the Digitest's autoclavable tip on the suspect tooth and press the "test" button. A gentle, pulsed stimulus begins to increase at a rate you've pre-selected (slow - medium - fast). Because there are no sudden jumps in power (as when you twist a knob), the patient feels the stimulus long before it crosses the "pain" threshold. Large LEDs indicate the strength of the stimulus.

FORAMATRON D-10 APEX LOCATOR

LOCATE THE APEX IN SECONDS THROUGH
ALL FLUID TYPES

\$895.00

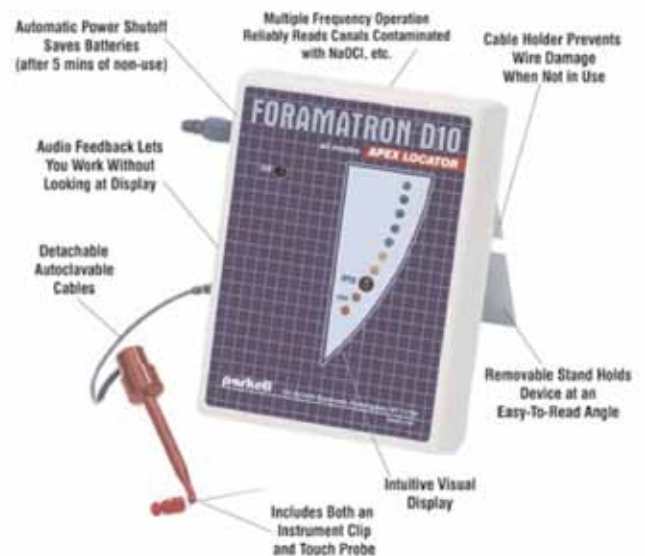

Turbo SENSOR Ultrasonic Scaler \$1,399.00



The TurboSENSOR's smart-circuitry senses whether you've inserted a 25KHz or 30KHz tip and automatically switches the operating frequency to match. It's done instantly, without your having to remember to flip a switch.

Higher "highs." Lower "lows." And greater power control in between. It gives you a higher maximum power setting for patients with extreme calculus. Yet for subgingival use, the tip begins vibrating at far lower power than any previous Parkell scaler. An orange LED tells you when you're operating in this low-power "perio" range. Turbo mode .. but without the Turbo button. To jump from low-power debridement to pedal-to-the-metal calculus-blasting, you simply increase foot-pressure, and the machine jumps into Turbo mode. Autoclavable, of course. Just unscrew the patented handpiece sheath and pop it into the autoclave. The TurboSENSOR comes with an extra sheath, so you can keep on scaling while one sheath is in the oven..

TurboSENSOR Includes two autoclavable handpiece sheaths, in-line water filter with spare filter disk, water connection with quick-connect. Accepts any 25KHz or 30KHz insert (Not included). Size: 1-3/4"H x 5-1/2"W x 7"D, Weight: 3 lbs.

SENSIMATIC ELECTROSURGERY 600SE


\$1,599.00



Solid state automatic circuitry for safety and performance. No vacuum tubes to break or burn out.

No warm-up time needed.

Size / weight: 3"H x 6-1/2"D x 11"W, 8.24 lbs. Power Output : 70 watts

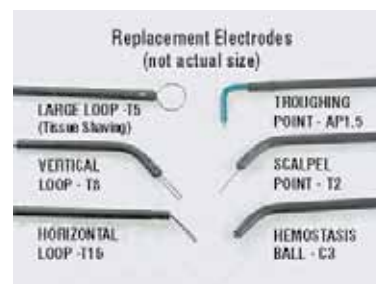

Supplied with 6 autoclavable electrode tips, foot controller, all accessories, and comprehensive operator's manual.

TIPS \$119.00 each 3 + \$109.00 ea